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
Product Introduction

EN

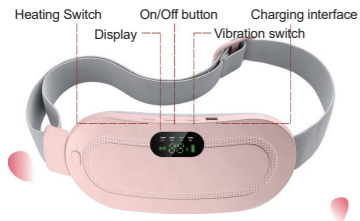Long press the power button to turn on or turn off the machine. Short press the warm switch to select temperature 122°F/131°F/140°F/STOP (50°C/55°C/60°C/STOP); Short press the vibration switch to select vibration mode- M1/M2 /M3/M4.

About this item

1.This warm heating pads uses the latest upgraded graphene technology. Fast heating portable cordless heating pad, it only takes 5 seconds to heat up and increase temperature. Helps relieve period cramp/abdominal/tummy/stomach or back pain.

2.The electric portable cordless heating pad has 4 vibration modes and 3 levels of heating. You can turn on at the same time. Allow you customize a setting for yourself to achieve the ultimate comfort. Three switchable temperatures, 122°F, 131°F and 140°F (50°C, 55°C and 60°C)

3.The portable heating pad is designed to automatically shut down. If you fall asleep during use, our product will automatically turn off within 30 minutes when it detects that you have not made any movement, which can ensure your safety. If you want to use it again, just turn it on again.

Specification parameters

Product Model	LJ-Y703
Input / Power	5V 1A / 3W
Battery Capacity	1800mAh
Charing Time	5 Hours

Charging/Recharging

- Please charge the heating pad before use.
- To charge the heating pad, plug the charger into an outlet and into the Charging interface.
- The indicator flashes red when charging, and will be solid green when complete.
- For safety considerations, the heating pad will NOT work while charging.
- For better maintenance, please recharge the heating pad once at least every 2 months if you don't use it regularly.

Safety Instructions


- Do NOT use while sleeping.
- Please avoid using the heating pad in the same position on your body for a long time and pay attention to your skin condition, if there is any abnormal on your skin, pls removed heating pad immediately.
- Do NOT bring your heating pad near fire or use all kind of corrosive liquid.

Warning

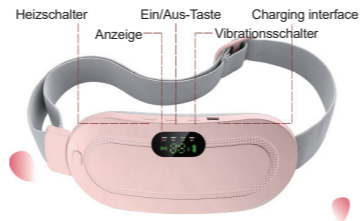- 1.Pregnant women, cervical tumor patients disabled.
- 2.Gout onset period or not completely healed with gout are forbidden to use.
- 3.Gravity disease and heart disease are forbidden to use.
- 4.The skin wound is still unhealed or patients with skin abnormalities are forbidden to use.

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