
### Warm reminder for product use

Dental experts suggest: the first time the user uses the gentle mode for one week, and then presses the switch button in the normal mode after getting used to it, it can be used, and the following points should be paid attention to when using:  
 (a) The direction of water flow is perpendicular to the teeth or gums  
 (b) Move slowly along the teeth and clean the teeth in order;  
 (c) The nozzle is aligned perpendicular to the gums;  
 (d) Thoroughly clean teeth, trays, crowns and dentures.  
 For the first use, do not spray water directly along the periodontal pocket for a long time, so as not to cause gum discomfort.

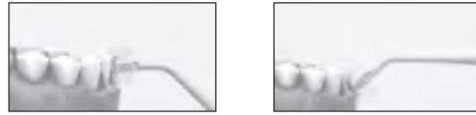
### Daily maintenance

A. If you need to clean the product, please wipe it with water or a damp cloth, and then use a dry cloth or other dry paper towels to wipe off the water left on the surface.  
 B. It is strictly forbidden to soak the product body in water.  
 C. Do not clean this product with hot water exceeding 50 degrees Celsius.  
 D. Do not bend, stretch or twist the nozzle hose.  
 E. It is recommended to replace the nozzle every 3 months.  
 F. If you do not use this product for a long time, you need to empty the water in the water tank and place it in a cool and dry place.


### Troubleshooting

Common problem	Possible Causes	Solution
Can't work or can't boot	New purchase or no more than three months	Please charge for more than 4 hours before use
The water pressure is very low during work	There is something blocked in the nozzle or the nozzle is deformed	Replace the nozzle
Can only be used for a short time after charging	The battery has reached the end of its useful life	Please contact customer service
The charging indicator does not light up or jumps randomly when charging	Insufficient charging time	Please charge for more than 4 hours before use
Do not spray water when working	The power is not connected	Please connect the power correctly
Common problem Can't work or can't boot	Adapter (charging head) does not match	Replace the adapter or plug directly into the USB socket
	Product circuit is faulty	Please contact customer service
The water pressure is very low during work	Newly purchased or left unused for a period of time	1. Pull out the nozzle to rinse the nozzle hole on the main body 2. Suck hard against the nozzle

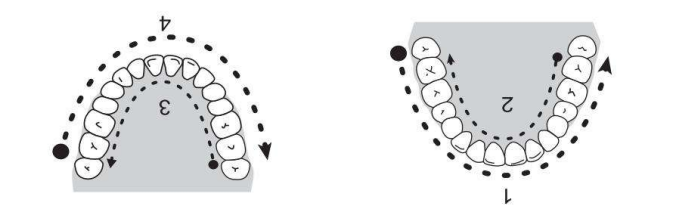flushing sequence for the other row of teeth is the same as above.  
 in the direction (such as 2 or 3) until it reaches the gum at the last tooth. The over-align the gums and continue to the same direction from the inside. Please as 1 or 4) until the gums of the last tooth are then turn the nozzle smoothly when flushing the teeth. After choosing the gear that suits you, align or the gums, and open the nozzle slightly so that the water flow can flow out the teeth, extend the nozzle into the mouth, aim at the gap between the teeth. Choose the gear that suits you. It is recommended to start the trial from the lowest gear and gradually find the gear that suits you best. When cleaning

### The right way to clean your teeth
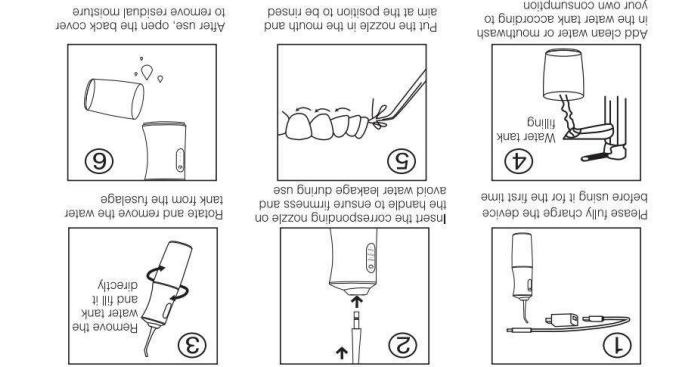
### Steps for usage

● Switch the mode to the pulse gear, and the indicator light also switches to the corresponding indicator light, and cycle in turn.  
 ● Short press to shut down.  
 ● When the battery of the instrument is charging, the three indicator lights will flash once, and when the battery is full, the three indicators will be on.  
 ● The instrument will automatically shut down after 5 minutes of operation.

### Environmental protection

Products used in this watch  
 Portable electric tooth flusher  
 Hazardous substance declaration  
 According to "Marking Requirements for Restricted Use of Hazardous Substances in Electrical and Electronic Products" GB/T 26572-2011 (also known as China RoHS2.0), the following section lists the names and contents of hazardous substances that may be contained in this product.

Part Name	Harmful materials					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated Diphenyl Ether (PBDE)
Plastic components	○	○	○	○	○	○
Hardware components	○	○	○	○	○	○
container	○	○	○	○	○	○
Motor components	×	○	○	○	○	○
battery	×	○	○	○	○	○
Data line	×	○	○	○	○	○
circuit board	×	○	○	○	○	○


This form is compiled according to the rules of SJ/T11364  
 ○: Indicates that the content of the hazardous substance in all homogeneous materials of the part is below the limit requirement specified in GB/T 26572-2011 (IEC62321).  
 ×: Indicates that the content of the hazardous substance in at least one of the homogeneous materials of the part exceeds the limit requirement of GB/T26572. However, the content of the parts marked with "x" in the above table exceeds the content because there is no mature alternative technology in the industry.

### Environmental protection use period

The mark on the product body indicates that the environmental protection use period is 10 years. The environmental protection use period of electrical and electronic products means that the toxic and hazardous substances or elements contained in electrical and electronic products will not leak to the outside or suddenly mutate, and users of electrical and electronic products will not be exposed to the environment when using the electrical and electronic products. The time limit for causing serious pollution or causing serious damage to the human body or property. During the environmental protection period, please use this product in accordance with the instruction manual. Remainder of the Regulations on the Recycling and Disposal of Waste Electrical and Electronic Products in order to better care and protect the earth, when users no longer need this product or the product life ends, please abide by the relevant national laws and regulations on the recycling and disposal of waste electrical and electronic products, and remove it. Hand it over to a local manufacturer with nationally recognized recycling and processing qualifications for recycling.

### Method of operation

- Short press to turn on, the default mode is strong, and the corresponding indicator lights up 2 seconds the switch mode, the mode switches to soft gear, and the indicator light also switches to the corresponding indicator light.
- Standard nozzle
- Periodontal pocket nozzle
- Special for periodontal pocket cleaning
- Orthodontic nozzles are suitable for cleaning braces
- Tongue scraping nozzle is dedicated to cleaning tongue coating
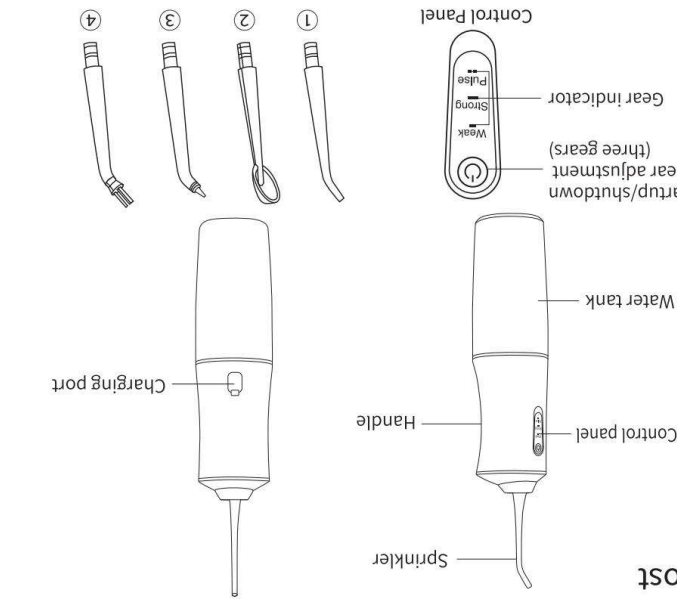
### How to dispose of rechargeable batteries

1. The built-in rechargeable lithium battery of this product contains substances that may pollute the environment. When disposing of this product, please do not put it together with general household garbage. The battery should be taken out.
2. Before disassembling this product and removing the lithium battery, be sure to turn off the power or cut off the power.
3. To remove the lithium battery, pry off the buckle on the top cover of the product, open the top cover of the product, and then take out the lithium battery.
4. The taken-out lithium battery should be safely disposed of and put in the battery recycling center designated by the government, which is conducive to environmental protection. Special Note: only take out the battery when discarding this product. Please make sure that the battery is exhausted before taking it out. Non-professionals should not disassemble it.

### Recycle

When disposing of recycled products, please don't throw them together with general household garbage, but take them to the designated recycling center. Everyone is responsible for protecting the environment!

#### Quality assurance

We hereby guarantee that the products are of high quality, reliable performance, safe and practical. In order to solve the customer's worries, our company implements responsible quality and after-sales service, which can be handled with this card.

1. If the main parts of the product fail during normal use within three months, repair and replacement parts are free of charge.
2. When repairing, you should bring the receipt of payment and warranty card to determine the warranty period.
3. The product warranty period is three months.

**Note**  
 1. Abnormal use and man-made damage are not covered by the warranty.  
 2. If the customer disassembles the machine by himself, the warranty qualification will be cancelled.  
 Date of purchase: \_\_\_\_\_  
 Where to buy: \_\_\_\_\_

#### Certificate of conformity

Inspection result: \_\_\_\_\_ Qualified \_\_\_\_\_

Inspectors: \_\_\_\_\_

Date of inspection: \_\_\_\_\_

Check  
03

## Portable electric tooth flusher

## user's manual

Disclaimer: Due to the iteration of product appearance, if there is a difference between the appearance of the pixels in this manual and the actual product, please refer to the actual product.

### Product parameter

- Product net weight: 355g
- Charging voltage: 5V
- Charging current: 1A
- Rated voltage: 3.7V
- Rated power: 3.5W
- Water tank volume: 230ml
- Battery capacity: 1200mAh

### Warning

- A. Do not use the product for other purposes other than oral cleaning.
- B. Do not use the nozzle to press the gums or teeth too hard to avoid damaging the gums.
- C. Do not immerse the product in water to avoid product failure.
- D. Do not use the product in a high temperature and humid environment.
- E. Do not continue to use the power cord if the power cord is damaged.

### Danger

- A. Please read this manual carefully before use to avoid injury or damage caused by improper operation.
- B. Do not disassemble this product privately to avoid product damage or electric shock.
- C. Do not use water above 40 degrees Celsius to fill the water tank to flush teeth to avoid burns.
- D. If the product falls into the water accidentally during charging, please unplug the power plug before removing the product, and wait until the water stains are completely dry before charging. E. If the power cord or plug is damaged, please replace the power adapter before use.
- E. Avoid placing the product in a high temperature place, do not throw the product into a fire, the battery may explode.

### Note

- A. Do not use this product without the nozzle installed.
- B. Some functional mouthwashes may be corrosive to certain extent. Please confirm the ingredients of the mouthwash before use, and do not use it indiscriminately.
- C. If children, the elderly or those who do not know how to use this product, please make sure that someone supervises and instructs the user.
- D. Quickly clean the oral cavity: Whether at home or eating out, it is the most effective way to quickly restore a clean oral cavity.
- E. It has a relatively good effect on the peripheral cleaning of orthodontic teeth, dental implants, and dental veneers.
- F. Help develop healthy daily oral hygiene habits and keep teeth fresh and clean at all times.

### Precautions

- A. Do not use strong gear to continuously wash damaged gums or bleeding teeth to prevent inflammation.
- B. Please adjust the flushing intensity of the pulse water flow in real time according to your own oral situation.
- C. Do not soak the whole machine in water, it may cause product failure or complete damage.
- D. If you do not use this product for a long time, please make sure that the water tank is dry and free of water before storing to prevent the growth of bacteria.

### Product utility

- A. The strong pulsed water flow can completely clean up bacteria, tooth stains, tartar and food residues in all parts of the oral cavity.
- B. The ultra-fine water stream can clean the gingival sulcus and deep between the teeth without barriers, and improve the health of the gums.
- C. It can effectively remove most of the subgingival plaque on the tooth surface, and the effect is much better than traditional dental floss.
- D. Quickly clean the oral cavity: Whether at home or eating out, it is the most effective way to quickly restore a clean oral cavity.
- E. It has a relatively good effect on the peripheral cleaning of orthodontic teeth, dental implants, and dental veneers.
- F. Help develop healthy daily interdental cleaning habits and keep teeth fresh and clean at all times.

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**Nota**  
 A. Do not use strong gear to continuously wash damaged gums or teeth.  
 B. If you do not use the instrument for a long time, please remove the battery from the instrument to avoid damage to the battery.

Consulte primero a un médico y confirme que es correcto antes de usarlo. Si los niños, los ancianos o aquellos que no saben cómo utilizar este producto, úsalo indiscriminadamente. Confirme los ingredientes del enjuague bucal antes de usarlo y no use el producto sin la boquilla instalada.  
 A. Do not use strong gear to continuously wash damaged gums or teeth.  
 B. If you do not use the instrument for a long time, please remove the battery from the instrument to avoid damage to the battery.

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**Peligro**  
 A. Do not use strong gear to continuously wash damaged gums or teeth.  
 B. If you do not use the instrument for a long time, please remove the battery from the instrument to avoid damage to the battery.

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**Produktparameter**

**Notiz**

- A. Do not soak the whole machine in water, it may cause product failure
- B. When using this product for a long time, please check the water tank regularly to ensure that the water level is sufficient and the water is clean.
- C. Please do not use the product for a long time continuously to prevent the motor from overheating.

**Notiz**

- A. Do not use strong gear to continuously wash damaged gums or teeth.
- B. Please adjust the timing intensity or the pulse water flow in real time according to your own oral situation.

**Notiz**

- A. Please read this manual carefully before use to avoid injury or damage to the product.
- B. Do not disassemble this product privately to avoid product damage or electrical shock.
- C. Please do not use this product in high temperature and high humidity environment to avoid product damage.
- D. Please do not use this product in water to avoid product failure.

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### Danger

### Précautions

frais et propre à tout moment.  
F. Help develop good oral hygiene habits and keep teeth healthy.  
E. Patients with periodontal disease or those who have had oral surgery in the past 2 months, or who have suffered from oral problems, please consult a doctor before use.  
D. Please keep it out of the reach of children and babies, and not use it as a toy.  
C. Si les enfants, les personnes âgées ou ceux qui ne savent pas comment utiliser ce produit, utilisez-le sans discernement.  
Veillez confirmer les ingrédients du bain de bouche fonctionnels peuvent être corrosifs dans une certaine mesure.  
B. Certains bains de bouche fonctionnels peuvent être corrosifs dans une certaine mesure.  
A. N'utilisez pas ce produit sans la buse installée.

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**Nota**  
 A. Do not soak the whole machine in water, it may cause product failure.  
 B. Please adjust the brushing intensity or the pulse water flow in real time in base alla propria situazione orale.  
 C. Do not use strong gear to continuously wash damaged gums or teeth, it may cause pain or discomfort.  
 D. Do not disassemble this product privately to avoid product damage or electro-shock.  
 E. Please read this manual carefully before use to avoid injury or damage to the device.

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