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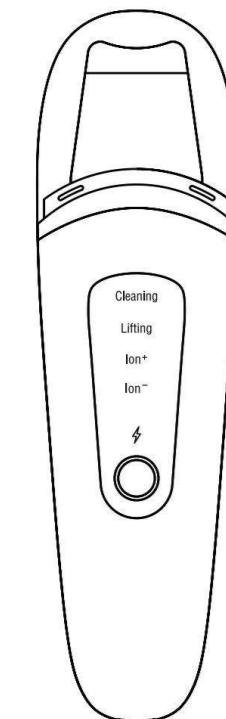
Ultrasonic Skin Scrubber

SS06

*Please read the instruction book carefully prior to use.

Safety precautions

Thank you for purchasing SS06 ultrasonic ion cleansing . Please read the instruction book carefully prior to use the product, and keep it properly so as to hand it over to any subsequent users.

About the instruction book: the product picture shown and related product feature introduction and the interpretation of some terms are only for reference. They are subject to the real product.

The instruction book contains safety instructions of proper use of the product. Therefore, please read the it carefully and keep it properly to avoid damage and danger caused by improper use. The product is personal can not be used for many people. All the safety instructions are for your safety

Important pre-use information

Such grade of mark indicates possible danger. It may cause personal injury or property loss in serious case. To this end, please comply with the corresponding warning description.

- If the charging line is damaged, please use the same specification of charging line to substitute it or contact after-sale department of our factory.
- If you feels ache or abnormal on the skin or your body, please stop using it immediately.
- When cleaning and maintaining the instrument, please close the machine and pull out the plug and do it in a shutdown state.
- Children are forbidden to use the product.
- Please don't use the product nearby flammable materials (alcohol, gasoline, nail polish remover, spray, etc.)
- Please don't immerse the product in the water or flush the product with tap water directly.
- If the product is damaged owing to careless drop or

01

collision, please don't use it again.

● Please don't dismantle or repair the product arbitrarily. Please contact the shop where you purchase the product for maintenance.

● Please use the charging line provided by our factory to ensure the voltage conform to the product specification requirements.

● The instrument must be cleaned in time after use to avoid accumulation of grease or other remnants.

● The products are forbidden to be used in the case of bath or shower, Please don't use the product nearby bathtub full of water, water basin or other vessels.

● Don't put the instrument in the dish-washing machine or microwave oven.

The following people are unsuitable to use the product


- Female in the pregnancy or menstruation period, or one with skin problem owing to hormone imbalance.
- Allergy population, person with skin disease, dermatitis, skin allergy and serious sunburn.
- Population with skin damage, eczema and swelling.
- Population during restoration after operation.
- Persons who have deficiencies in limb, sense, or spirit, or lack experience is forbidden to use the product unless under the supervision and guidance of guardian and the guardian is liable for the safety of user.

Feature

1. Help remove oils and other dirt.
2. Help loosen impacted oils and blackheads.
3. Help deliver your serums and moisturizers, optimizing the performance of skincare products.
4. Offer a comfortable experience on the face from a warm massage.

02

Name of parts

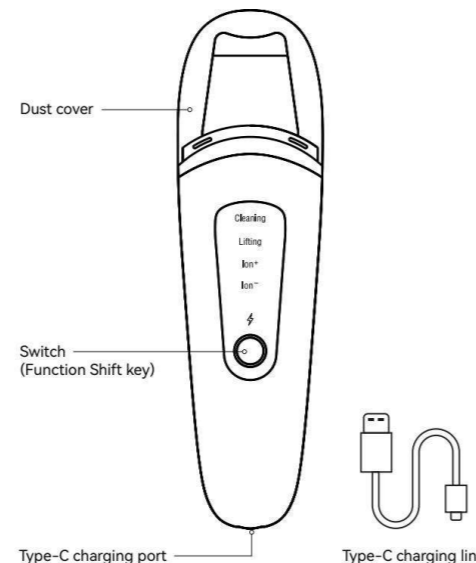Use instruction

1. Long press the button for 2 seconds to power on
2. Short press the button to switch gears
3. The red light flashes when charging, and the blue light is fully charged

03


Warm tips

1. Please charge it for more than 5hr before the first use
2. Please wet your face before use.If the skin is not thoroughly cleaned, it will not be able to absorb the essence well.
- 3.The recommended use frequency is 2-3 times per week to avoid overuse.
4. Allergic skin should reduce the time and frequency of use.
5. Do not rinse it directly under the tap, because it may damage the machine, please use a wet soft mop with soap and water to clean the stainless steel spatula
6. Please fully charged it and store the device in a clean, dry, and cool place.
- 7.The skincare cycle is 28 days, please use it 2 or 3 times a week to achieve better effects.

Steps for usage


Cleanse the face

Step 1:
Clean your face with facial foam, use a steamer/ warm towel / spray to keep your face wet.

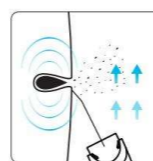Enjoy the massage


Step 2:
Using essence, oil, or serum and enjoy a comfortable experience on the face from the warm massage.

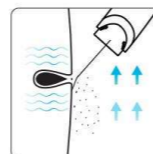04


Deeply cleanse pores

Step 3:
Keep your skin wet, Say goodbye to blackheads and clogged pores with this skin-cleansing tool.

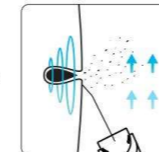Optimize skin care products

Step 4:
Using essence, lotion, or face mask, optimizes the performance of your skincare products.

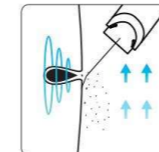Cleaning and maintenance

- Please keep the product away from humidity, and don't put the instrument in the dish-washing machine or microwave oven. Prior to cleaning, please close the machine.
- Please don't use sharp, corrosive, or volatile diluent.
- The product can be wiped gently with cold water or warm water or soft cloth or cotton pad.
- Please don't flush the instrument with water directly, because it may damage the machine.
- If the machine keeps idle for a long time (about three months), please fully charge the machine and then store the machine for fear of machine damage owing to battery power deficiency.

05

Troubleshooting


Fault	Possible reason	Measures
Shovel doesn't vibrate.	Shovel deformation.	Hold the shovel with both hands and correct it.
	Internal wire falls off.	Please contact after sale department of our factory.
Fail to start	Reduce the operation force, just touch the shovel gently.	Reduce the operation force, just touch the shovel.
	Battery power shortage.	Charge the machine for 2h prior to power on.
	Battery damage caused by too long resting period of the machine.	Please contact after-sale department of our factory.
Ion function failure	Failure caused by water entering the key.	Please contact after sale department of our factory.
	Hands don't touch the bus bar on the lateral side of machine.	During operation, the hands should touch the bus bar on the lateral side of machine.
	Internal circuit of machine is damaged owing to water ingress.	Please contact after-sale department of our factory.

*If the product cannot work normally, please stop using it and contact after-sale department of the factory.

Technical parameter

Working voltage: 4.2~3.3
Working current: 500mA
Blade Frequency: 24kHz~27kHz
Battery capacity: 500mah
Product Size: 173*48*25.5mm
Package Size: 187.5*69.5*33.5mm

06

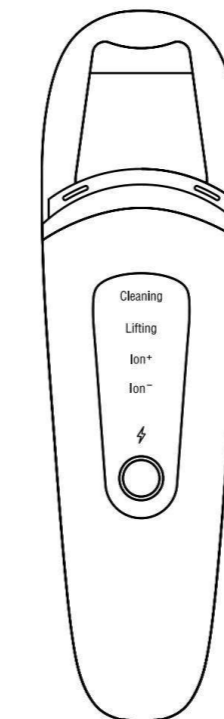Ultrasonic Skin Scrubber

SS06

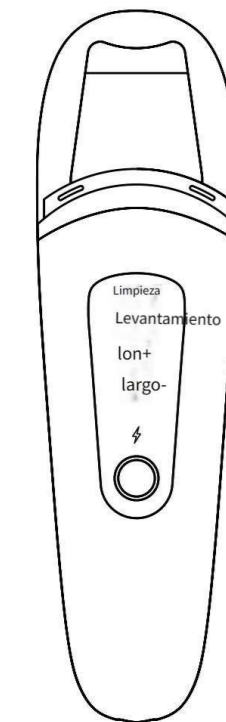
*Please read the instruction book carefully prior to use.

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