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
Made in China

Neck Beauty Device

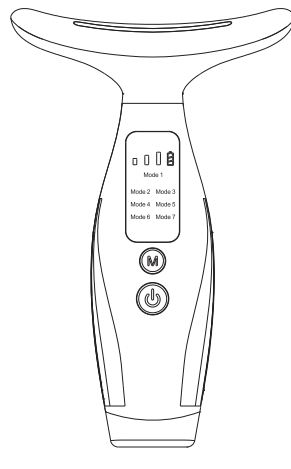Product Manual
Thank you very much for purchasing and using this product

Product Specifications


Product name	Neck beauty device
Product number	HY69
Product Size	139*90*38MM
Rated power	3W
Power adapter	DC5V \equiv 1A(USB)
Charging time	3 hours

Button definition

Message area
Long press the power on/off button: power on/off
Short press the mode key
Red light + hot compress + vibration, display mode 1
Orange light + hot compress + vibration, display mode 2
Green light + hot compress + vibration, display mode 3
Yellow light + hot compress + vibration, display mode 4
Blue light + hot compress + vibration, display mode 5
Blue light + hot compress + vibration, display mode 6
Purple light + hot compress + vibration, display mode 7

Power button
Seven modes switching work
Automatically shuts down after 15 minutes without operation

Motor strength
Short press the power button: the motor strength can be switched between levels 1-3.
The corresponding gear light is on
Plug in the charging cable when the phone is off or on
Charging status, battery indicator light flashes, battery is full
The indicator light is always on

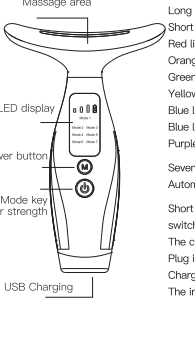
Product Maintenance

1. Please try to keep this product in a dry and room temperature environment or in a storage box.
2. Do not overcharge to avoid damaging the battery life.
3. Please wipe the product clean after charging and put it in the storage box if not used for a long time.
4. The outside of the product can be wiped clean with alcohol or diluted detergent.
5. During use, the device will become slightly hot, but it is normal for it not to be hot.
6. It is forbidden to rinse the fuselage directly with water.
7. If it is not used for a long time, please charge the battery at least once every 3 months.

Usage warning

1. It is prohibited to use it with medical electrical appliances implanted in the body, such as pacemakers, as it may cause equipment failure or trigger danger.
2. Do not use this product if your palms have skin problems such as wounds, eczema, blisters, and peeling.
3. Do not use this product in places with high humidity such as bathrooms.
4. Please do not use this product for other purposes than its intended use.
5. Please do not disassemble, repair or modify this product without permission. May cause equipment malfunction or cause injury.
6. Please do not give this product to children under six years old to avoid personal injury.

Packing list

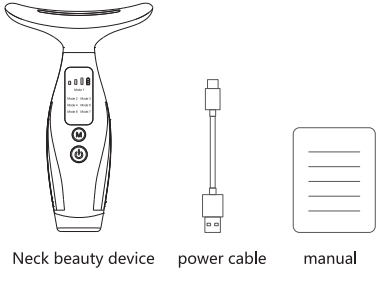Neck beauty device power cable manual

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