PARAMETER

Weighing range: 0.2kg–180kg

D=50g unit conversion: kg/lb

Power supply: 3V battery: 2 AAA batteries

Power off current: $\leq 25\text{mA}$ standby current: $\leq 10\text{mA}$

Low battery indicator light: [Lo] battery needs to be replaced

Overload indication: [o-Ld] Overload exceeds 183kg

Best operating temperature: 5–35°C

OKOK App Weighing instructions

Connect the APP for weighing, generate a report on the measurement results, and provide health recommendations. The data will be stored in the historical records. The trend of weight changes, including weekly, monthly, and quarterly reports, is represented by curves to show the changes in weight.

Measurement items: weight, body shape, obesity BMI, Multiple data including muscle mass, water content, basal metabolism, bone mass, protein, etc

User registration, login, password change, password retrieval and other functions

Manage user information, including nickname, birthday, height, gender Wait, add family members

precautions


1. Please keep the scale away from damp ground. If the glass is wet, do not use the scale to prevent slipping.
2. Please stand in the middle of the scale to prevent tipping over.
3. The scale panel is made of glass, please avoid impact.
4. Keep the surrounding environment dry and do not immerse it in water.
5. When the battery is in a low battery state, please remove the battery to prevent leakage and corrosion of the scale.

instructions

Step 1: Install the battery:

1. Open the battery cover
2. Install the battery and pay attention to the positive and negative poles
3. Close the lid
4. Press the [unit] button on the back to convert kg/lb

Tip: If not used for a long time, please remove the battery

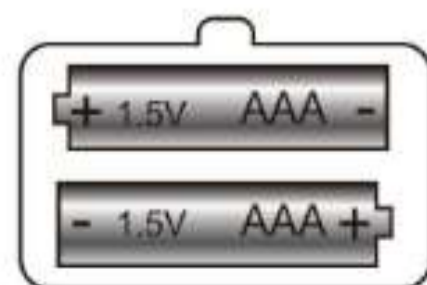Step 2: Install the OKOK application

1. Scan the QR code to download the app or search for the 'OKOK' app name in the app store to download and install it
2. Enable the BT, location, and nearby device permissions on the phone.
iPhone: Settings – Privacy – OKOK – Enable app permissions
Android: Settings – Privacy/Application Permissions – OKOK – Enable Permissions
3. The mobile network is normal

OKOK App

IOS7.0

Android4.3

Step 3: Register an account and connect to the APP

1. Open the app to start pairing devices.
2. Register as a user.
3. Fill in personal information correctly


Step 4: Measure and Weigh

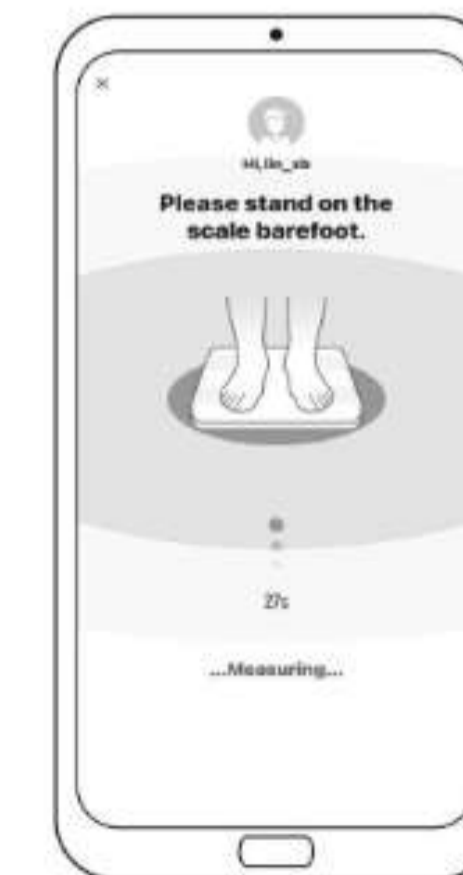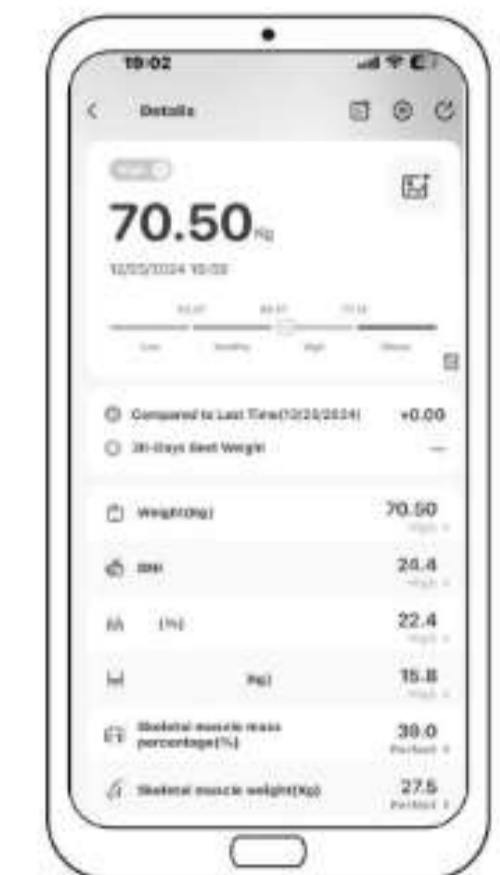
1. Place the scale on a hard and level surface, press and turn it back to zero.
2. Enter the "Dynamic Weighing" interface of the app
3. Stand on the scale and synchronize your weight to upload to the OKOK App

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